


Purine alkaloid

No.03217

- 1. Theobromine
- 2. Theophylline

- 3. Caffeine
- 4. Phenol

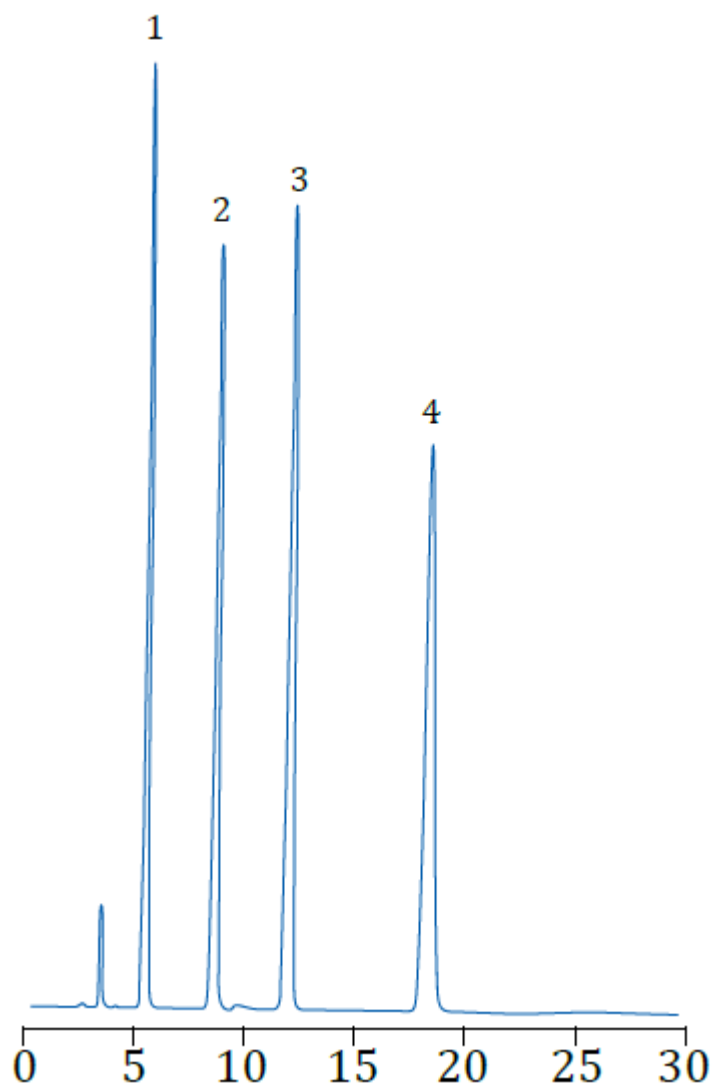

Column: Athena C18 4.6 × 150mm, 5µm (LAEQ-461571)

Mobile phase: methanol / water (25/75)

Flow rate: 1.0 mL/min

Detection: 254 nm

Column temperature: 40 °C

ATHENA C18 APPLICATION & ORDERING INFORMATION

Ordering Information:

Particle size diameter × length Item No.

Packings	Particle size	diameter × length	Cat. No
Athena C18	5µm	2.1 × 50mm	8.210571.0001
Athena C18	5µm	2.1 × 100mm	8.211071.0001
Athena C18	5µm	2.1 × 150mm	8.211571.0001
Athena C18	5µm	2.1 × 200mm	8.212071.0001
Athena C18	5µm	2.1 × 250mm	8.212571.0001
Athena C18	5µm	4.6 × 50mm	8.460571.0001
Athena C18	5µm	4.6 × 100mm	8.461071.0001
Athena C18	5µm	4.6 × 150mm	8.461571.0001
Athena C18	5µm	4.6 × 200mm	8.462071.0001
Athena C18	5µm	4.6 × 250mm	8.462571.0001
Athena C18	10µm	4.6 × 150mm	8.461574.0001
Athena C18	10µm	4.6 × 250mm	8.462574.0001

ATHENA C18 APPLICATION & ORDERING INFORMATION

www.doalab.com

info@doalab.com (+90) 212 509 21 12