


ATHENA C18-WP APPLICATION & ORDERING INFORMATION

Procainamide

LC-10039

1.	Uracil	66-22-8
2.	Procainamide	51-06-9
3.	N-Acetylprocainamide	32795-44-1

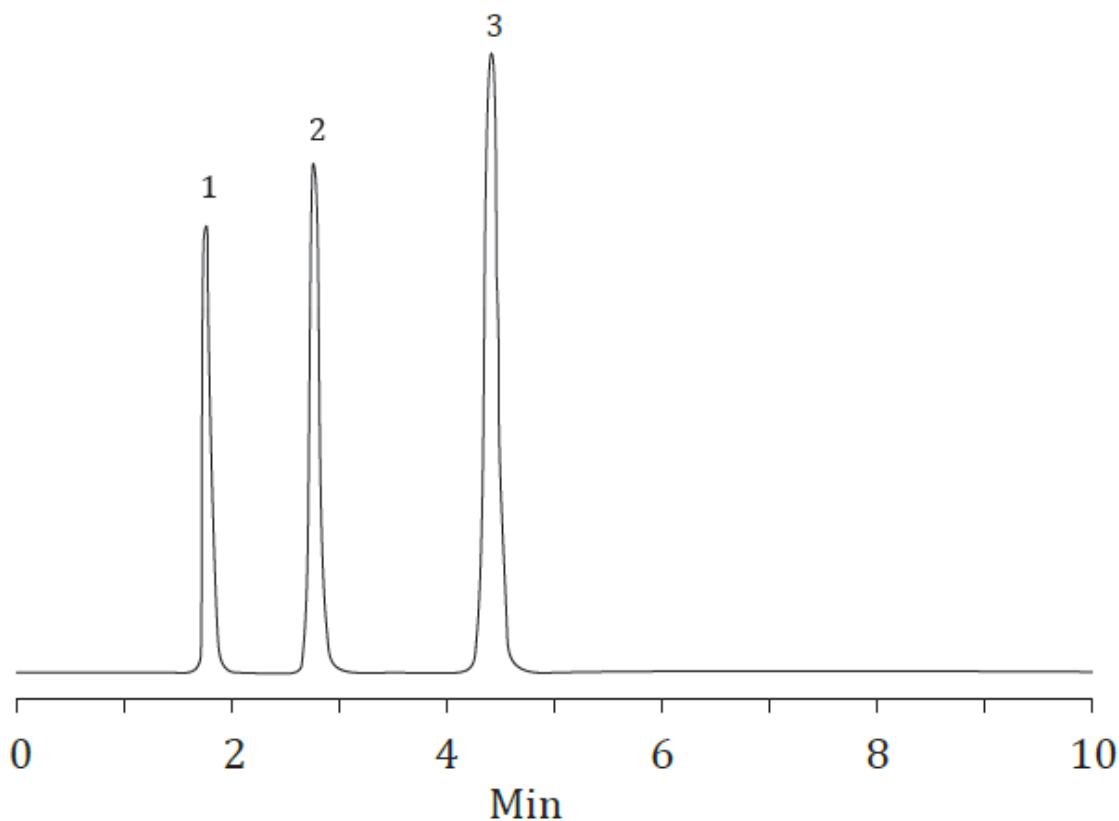

Column	Athena C18-WP 4.6 × 150mm, 5 μ m (LAEQ-461572)
Mobile Phase	Methanol /20mM KH ₂ PO ₄ - K ₂ HPO ₄ (pH 7.0) (40/60)
Flow Rate	1.0 ml/min
Column Temperature	40 °C
Detector	254 nm

Ordering Information:

Packings	Particle size	diameter × length	Cat. No
Athena C18-WP	3µm	2.1 × 50mm	8.210573.0001
Athena C18-WP	3µm	2.1 × 100mm	8.211073.0001
Athena C18-WP	3µm	2.1 × 150mm	8.211573.0001
Athena C18-WP	3µm	2.1 × 200mm	8.212073.0001
Athena C18-WP	3µm	2.1 × 250mm	8.212573.0001
Athena C18-WP	3µm	4.6 × 50mm	8.460573.0001
Athena C18-WP	3µm	4.6 × 100mm	8.461073.0001
Athena C18-WP	3µm	4.6 × 150mm	8.461573.0001
Athena C18-WP	3µm	4.6 × 200mm	8.462073.0001
Athena C18-WP	3µm	4.6 × 250mm	8.462573.0001
Athena C18-WP	5µm	2.1 × 50mm	8.210572.0001
Athena C18-WP	5µm	2.1 × 100mm	8.211072.0001
Athena C18-WP	5µm	2.1 × 150mm	8.211572.0001
Athena C18-WP	5µm	2.1 × 200mm	8.212072.0001
Athena C18-WP	5µm	2.1 × 250mm	8.212572.0001
Athena C18-WP	5µm	4.6 × 50mm	8.460572.0001
Athena C18-WP	5µm	4.6 × 100mm	8.461072.0001
Athena C18-WP	5µm	4.6 × 150mm	8.461572.0001
Athena C18-WP	5µm	4.6 × 200mm	8.462072.0001
Athena C18-WP	5µm	4.6 × 250mm	8.462572.0001

ATHENA C18-WP APPLICATION & ORDERING INFORMATION

www.doalab.com

info@doalab.com (+90) 212 509 21 12