


Pigment**LC-20080**

1.	Amaranth	915-67-3
2.	Ponceau 4RC	2611-82-7
3.	Sunset Yellow	2783-94-0
4.	Brilliant Blue	3844-45-9
5.	Erythrosin B disodium salt	16423-68-0

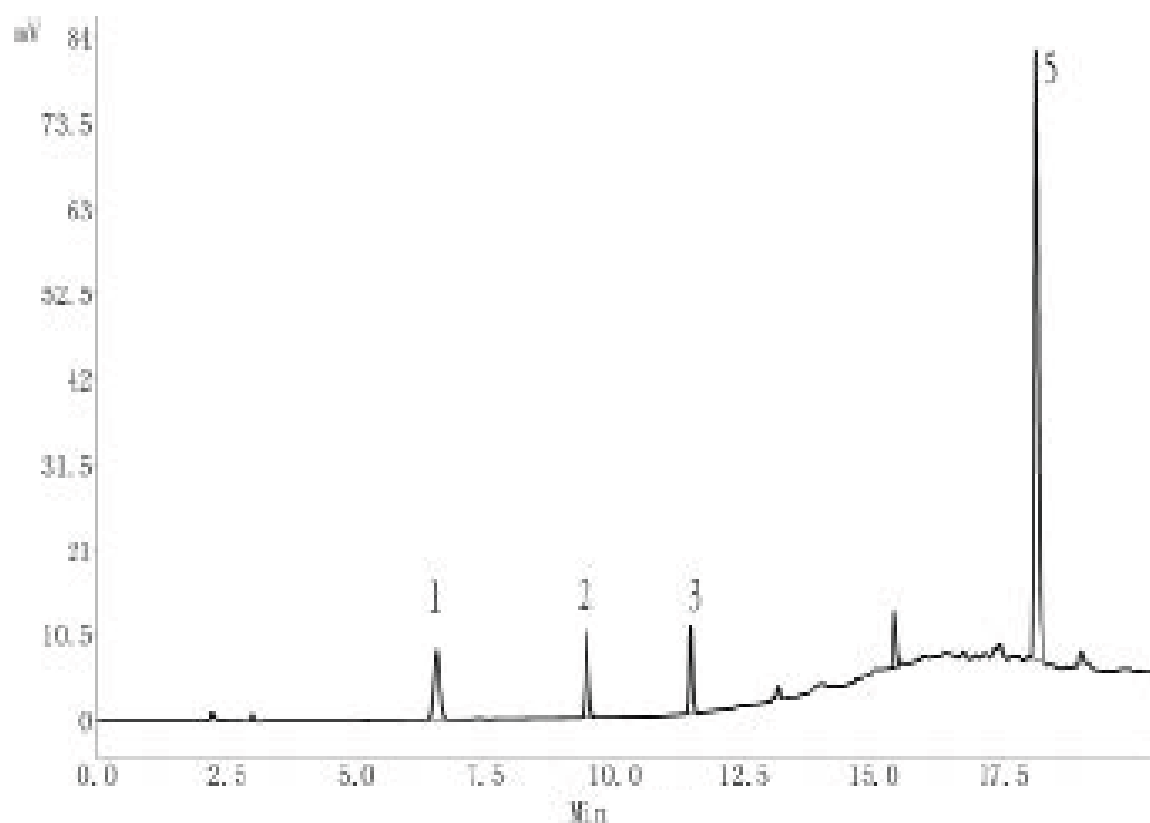

Column	Athena C18-WP 4.6mm*250mm, 5um (LAEQ-462572)		
Mobile Phase	A; Methanol B; Ethanol amine		
	T (min)	Methanol	Ammonium acetate
	0	20	80
	5	35	65
	12	98	2
	18	98	2
Flow Rate	1.0ml/min		
Column Temperature	30°C		
Detector	254nm		