


Tocopherol isomers

No.03227

1. α -Tocopherol
2. β -Tocopherol
3. γ -Tocopherol
4. δ -Tocopherol

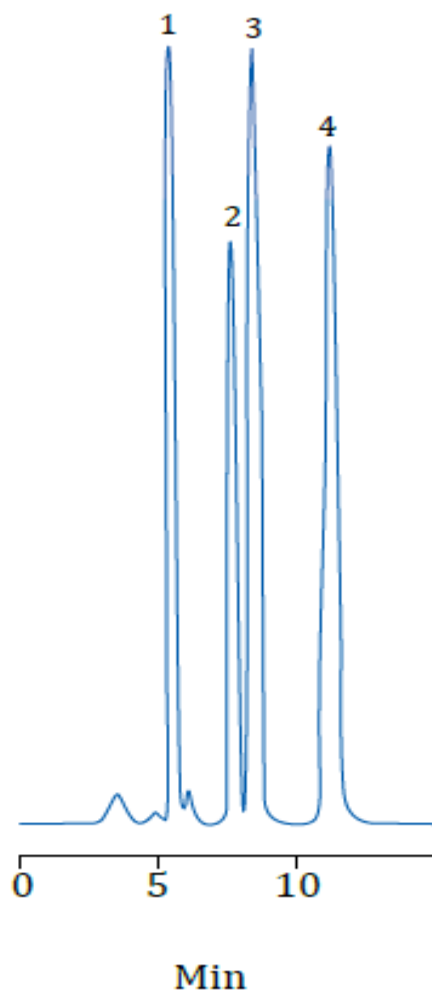

Column: Athena NH₂ 4.6 x 150mm, 5 μ m (8.461576.0001)

Mobile phase: n-hexane / ethyl acetate (70/30)

Flow rate: 1.0 mL/min

Detection: 295 nm

Column temperature: 40 °C

temperature:

ATHENA NH APPLICATION & ORDERING INFORMATIONS

Ordering Information:

| Packings | Particle size | diameter × length | Cat. No |
|------------|---------------|-------------------|---------------|
| Athena NH2 | 3µm | 2.1 × 50mm | 8.210567.0001 |
| Athena NH2 | 3µm | 2.1 × 100mm | 8.211067.0001 |
| Athena NH2 | 3µm | 2.1 × 150mm | 8.211567.0001 |
| Athena NH2 | 3µm | 2.1 × 200mm | 8.212067.0001 |
| Athena NH2 | 3µm | 2.1 × 250mm | 8.212567.0001 |
| Athena NH2 | 3µm | 4.6 × 50mm | 8.460567.0001 |
| Athena NH2 | 3µm | 4.6 × 100mm | 8.461067.0001 |
| Athena NH2 | 3µm | 4.6 × 150mm | 8.461567.0001 |
| Athena NH2 | 3µm | 4.6 × 200mm | 8.462067.0001 |
| Athena NH2 | 3µm | 4.6 × 250mm | 8.462567.0001 |
| Athena NH2 | 5µm | 2.1 × 50mm | 8.210577.0001 |
| Athena NH2 | 5µm | 2.1 × 100mm | 8.211077.0001 |
| Athena NH2 | 5µm | 2.1 × 150mm | 8.211577.0001 |
| Athena NH2 | 5µm | 2.1 × 200mm | 8.212077.0001 |
| Athena NH2 | 5µm | 2.1 × 250mm | 8.212577.0001 |
| Athena NH2 | 5µm | 4.6 × 50mm | 8.460577.0001 |
| Athena NH2 | 5µm | 4.6 × 100mm | 8.461077.0001 |
| Athena NH2 | 5µm | 4.6 × 150mm | 8.461577.0001 |
| Athena NH2 | 5µm | 4.6 × 200mm | 8.462077.0001 |
| Athena NH2 | 5µm | 4.6 × 250mm | 8.462577.0001 |

ATHENA NH APPLICATION & ORDERING INFORMATIONS

www.doalab.com

info@doalab.com (+90) 212 509 21 12